


Hydraulic Surgical Table

Product introduction

YA-04M is the mobile surgical table is fully manual operated for easy and effective operation in order to facilitate access for either surgical team or C-arm and ensure patient safety and comfort.

The surgical tables are hydraulically operated in elevation and crank operated in all of tabletop articulations through head-end control position.

This avoids wasting time looking for controls located over the table, or distributing the surgical team when the table needs to be adjusted during any surgical procedures.

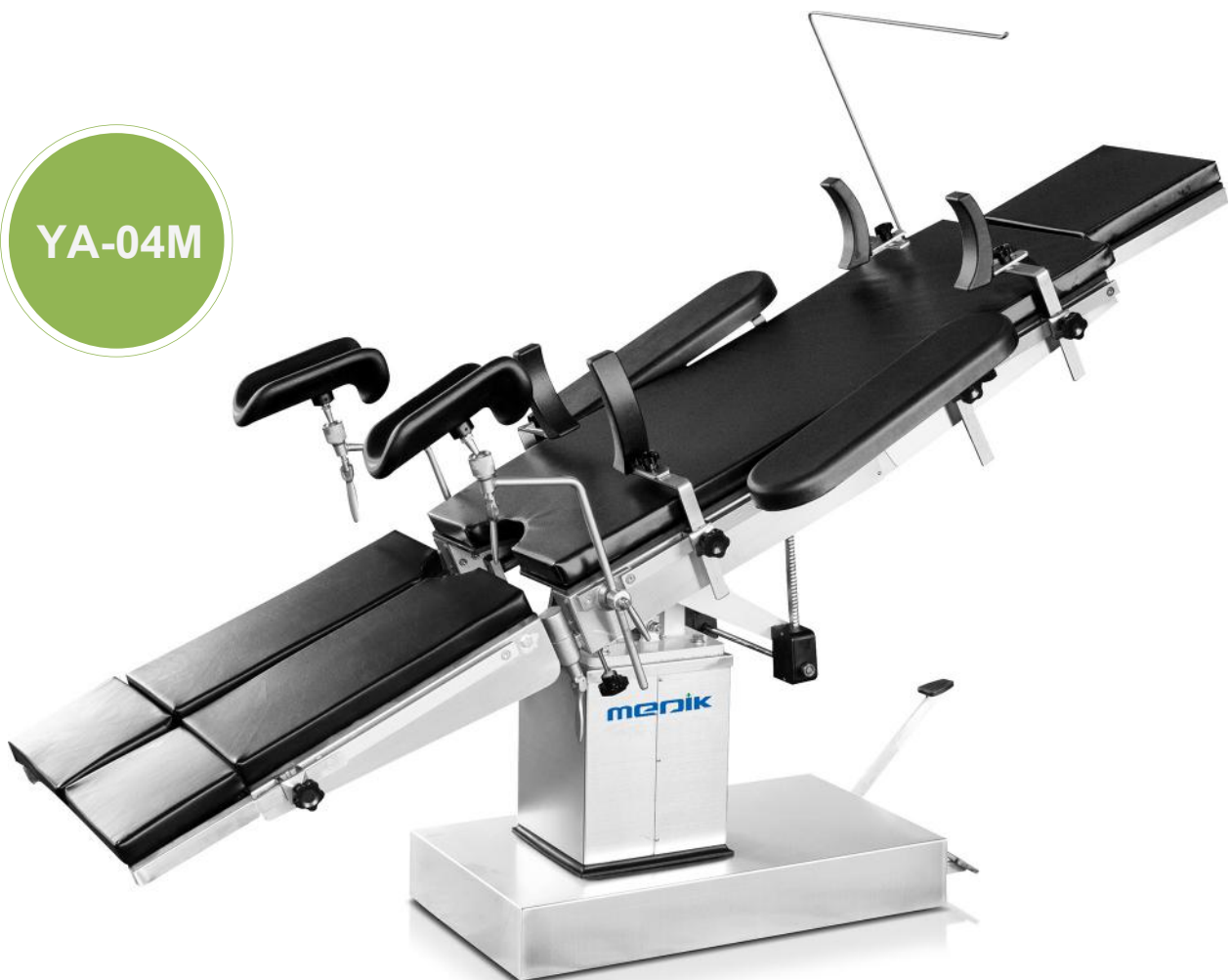

Hydraulic Surgical Table

Main details



Main material

The base material is stainless steel304, which is easier to clean, modern and stable.



Manual control system

various actions are controlled by hydraulic pneumatic system and mechanical driven system.



Kidney bridg

built-in kidney bridge, can be adjusted up to 100mm through the crank.




Brake system

The table can be moved easily or in a stable position thanks to the special foot brake system.



Height adjust

Footstep hydraulic elevating, imported Y type sealing ring and is durable

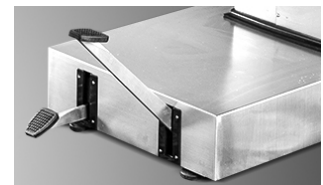

Technical description

and Table specifications

Table specifications

Reverse trendelenburg

0-25°



Trendelenburg

0-20°



Lateral tilt right

0-18°



Lateral tilt left

0-18°



Lowest position of table top

750 mm



Highest position of table top

1000 mm



Back section up motion

0-60°



*The page shows the main parameters, Please contact with us for full parameters or more info.

Technical description

and Table specifications


Standard accessories



Body support



Leg support

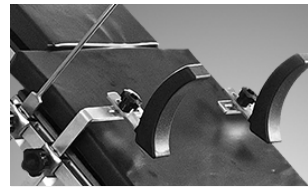

Shoulder support



Arm board

Technical parameter

	PARAMETER
Overall size	2100*500*750-1000
Reverse trendelenburg	mm
Trendelenburg	25°
Lateral tilt left/right	20°
Head section up/down	18°/18°
Back section up/down	30°/90°
Leg section down	60°/20°
Kidney elevator	90°
mechanism	100mm
Slide motion	300mm
Loading capacity	200kg

Technical configuration

	PARAMETER
Bed frame	1piece
Mattress	1set
Anesthesia screen	1piece
Shoulder supportr	1pair
Arm board	1pair
Body support	1pair
Leg support	1pair
Foot plate	1pair
Kidney bridge handle	1piece
Fixing clamp	8pieces
Oil pot	1piece

*The page shows the main parameters, Please contact with us for full parameters or more info.

ABOUT US

Medik is a group focusing on surgical room solution and nursing care solution.

We started to manufacture surgical table&lamp since 2011 with a professional team consisting of experienced designer and engineer from international top brand. With 10 years accumulation, Medik is becoming a new star to provide reliable product and valuable solution for surgical room.



YOU TUBE



MEDIK

Pls scan the QR code for more info

info@medikmedical.com | +86 13962231292

Medik reserves the right to make changes without notice in design, specifications and models. The only warranty Medik makes is the express written warranty extended on the sale or rental of its products.

© 2022 Medik Services, Inc. ALL RIGHTS RESERVED. APR53703 rev 1 01-JAN-2022 ENG – US