

# Electric Operating Table

## Product introduction

Engineered for increasingly complex surgical scenes, YA-05E comes with intelligent automatic control system satisfying various demands of general surgeries.

The structure enables it to be adapted to match the size of the patient and in no time at all it can be equipped for any surgical discipline. 350mm long horizontal sliding that can be suitable for both X-ray and C-arm.

YA-05E is an ideal electric operating table.




# Electric Operating Table

## Main details



### Main material

The base material is stainless steel304, which is easier to clean, modern and stable.

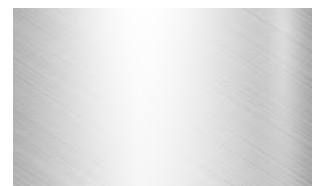

### Handle controller and column controller

Equipped with hand controller and column controller, the icons are clear and easy to operate.



### Kidney bridge

built-in kidney bridge, can be adjusted up to 110mm through the crank.



### 5" Castors central brake

Central control casters, large wheels with adiameter of 125mm, beautiful and easy to move.



### Power

Power with the fuse and with a master switch, it needs to be pressed when using internal battery maintaining power supply for 50-80 operations.



# Technical description

## and Table specifications

### Table specifications

Lowest position of table top

**710** mm



Highest position of table top

**1010** mm



Back section up motion

**0-80°**



Back section down motion

**0-20°**



Reverse trendelenburg

**0-25°**



Trendelenburg

**0-25°**



Lateral tilt right

**0-15°**



Lateral tilt left

**0-15°**



Slide motion

**350** mm



Head section up

**0-40°**



Head section down

**0-90°**



Leg Section Rest

**-90-15°**



\*The page shows the main parameters, Please contact with us for full parameters or more info.

# Technical description

## and Table specifications

### Standard accessories



**Body support**



**Leg support**



**Arm board**



**Shoulder support**

### Technical parameter

	PARAMETER
Overall size	2100*550*710-1010mm
Reverse trendelenburg	25°
Trendelenburg	25°
Lateral tilt left/right	15°/15°
Head section up/down	45°/90°
Back section up/down	80°/20°
Leg section up/down	15°/90°
Leg section outwards	90°
Kidney elevator mechanism	110mm
Slide motion	350mm
Loading capacity	250kg

### Technical configuration

	PARAMETER
Bed frame	1piece
Mattress	1set
Anesthesia screen	1piece
Shoulder supportr	1pair
Arm board	1pair
Body support	1pair
Leg support	1pair
Foot plate	1pair
Hand remote	1piece
Fixing clamp	8pieces

\*The page shows the main parameters, Please contact with us for full parameters or more info.

## ABOUT US

Medik is a group focusing on surgical room solution and nursing care solution.

We started to manufacture surgical table&lamp since 2011 with a professional team consisting of experienced designer and engineer from international top brand. With 10 years accumulation, Medik is becoming a new star to provide reliable product and valuable solution for surgical room.



YOU TUBE



MEDIK

**Pls scan the QR code for more info**

**[info@medikmedical.com](mailto:info@medikmedical.com) | +86 13962231292**

Medik reserves the right to make changes without notice in design, specifications and models. The only warranty Medik makes is the express written warranty extended on the sale or rental of its products.

© 2022 Medik Services, Inc. ALL RIGHTS RESERVED. APR53703 rev 1 01-JAN-2022 ENG – US