

Hydraulic Surgical Table

Product introduction

YA-06M is a hydraulic and mechanical operated mobile operating table with simple and effective operation to ensure patient safety and comfort.

The operating table is adjusted by ratchet operation on both sides, and the hydraulic pump can be raised and lowered by foot control, which can avoid wasting time to find the controller located above the operating table.



Hydraulic Surgical Table

Main details



Main material

The base material is stainless steel304, which is easier to clean, modern and stable.



Manual control system

Various actions are controlled by hydraulic pneumatic system and ratchet system.



Brake system / Height adjust

The table can be moved easily or in a stable position thanks to the special foot brake system. Footstep hydraulic elevating, imported Y type sealing ring and is durable.



Lateral tilt left / right

The crank on the right controls the left and right tilt of the table top.



Kidney bridg


Built-in kidney bridge, can be adjusted up to 100mm through the crank.



Technical description

and Table specifications

Table specifications

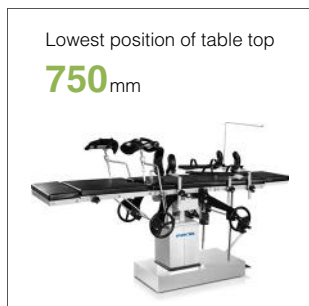

*The page shows the main parameters, Please contact with us for full parameters or more info.

Technical description

and Table specifications

Standard accessories



Shoulder support



Leg support



Arm board



Body support

Technical parameter

	PARAMETER
Overall size	2050*500*750-1000
Reverse trendelenburg	mm30°
Trendelenburg	20°
Lateral tilt left/right	18°/18°
Head section up/down	45°/90°
Back section up/down	75°/18°
Leg section down	90°
Kidney elevator mechanism	110mm
Loading capacity	250kg

Technical configuration

	PARAMETER
Bed frame	1piece
Mattress	1set
Anesthesia screen	1piece
Shoulder supportr	1pair
Arm board	1pair
Body support	1pair
Leg support	1pair
Foot plate	1pair
Kidney bridge handle	1piece
Fixing clamp	8pieces
Oil pot	1piece

*The page shows the main parameters, Please contact with us for full parameters or more info.

ABOUT US

Medik is a group focusing on surgical room solution and nursing care solution.

We started to manufacture surgical table&lamp since 2011 with a professional team consisting of experienced designer and engineer from international top brand. With 10 years accumulation, Medik is becoming a new star to provide reliable product and valuable solution for surgical room.



YOU TUBE



MEDIK

Pls scan the QR code for more info

info@medikmedical.com | +86 13962231292

Medik reserves the right to make changes without notice in design, specifications and models. The only warranty Medik makes is the express written warranty extended on the sale or rental of its products.

© 2022 Medik Services, Inc. ALL RIGHTS RESERVED. APR53703 rev 1 01-JAN-2022 ENG – US