


MK-E370JL

Mobile Examination Light

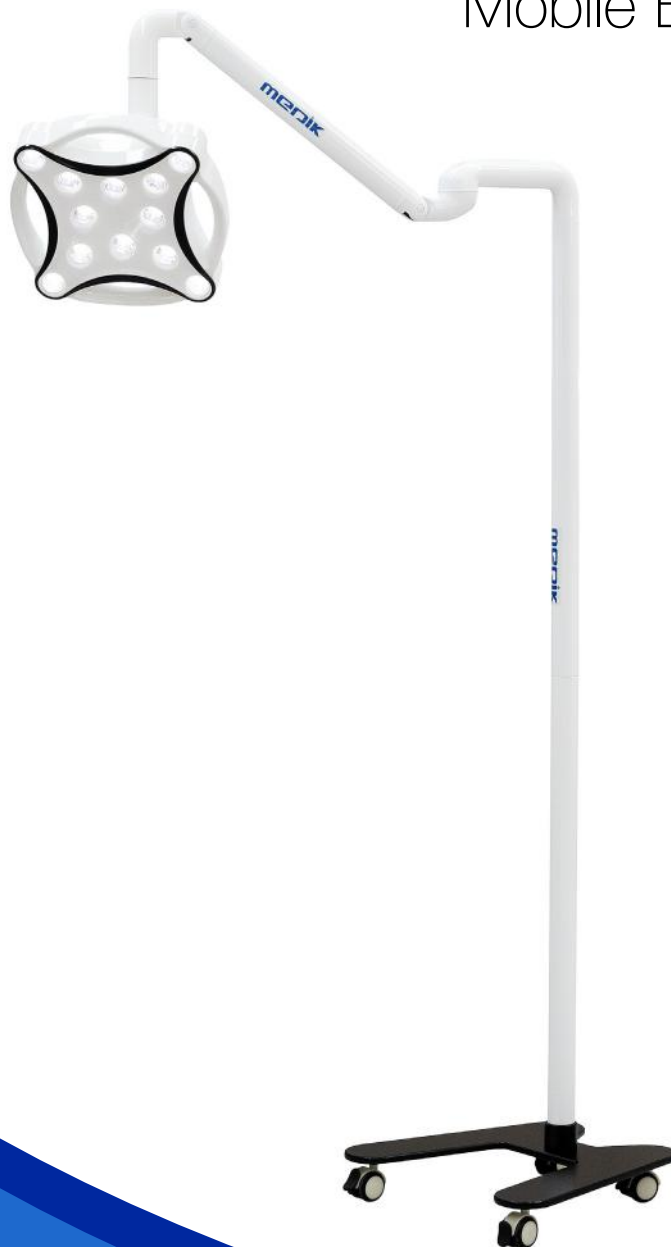

*Better care
Better comfort!*

MK-E370JL

MOBILE EXAMINATION LIGHT

DESCRIPTION

Flexible hose with two joints to get wider moving space for the light head, elegant appearance and excellent light parameters, all the features make the MK-E370JL Mobile Examination Light is a perfect solution designed for patient exam, minor surgeries in the hospitals.

FEATURES

- 1.Gesture-Controlled Brightness Adjustment
- 2.User-Friendly Adjust Brightness Button
- 3.Designed for hospital mobility, this examination light offers easy maneuverability to different locations.
- 4.Designed for enhanced surgical assistance

STANDARD SPECIFICATION

Power	30W
Luminance (LUX)	5200-100000lux
Spot diameter(mm)	130mm
Color temperature (K)	4000-5000
LED life	50000hrs



LED Bulb

Powered by long-lasting LED bulbs, ensures an extended lifespan and minimal heat emission. Experience efficiency and durability.



Sensing adjustment brightness

Effortlessly control the brightness of the light using hand gestures. The intuitive gesture-sensing technology allows doctors to adjust the illumination without physical contact, thereby reducing cross-infection.



Multi-segmented arms

Featuring multi-segmented arms allows flexible adjustment at various angles. Tailor the illumination precisely to your needs, offering adaptability for diverse applications.



Clear switch

Equipped with a master switch and brightness adjustment buttons on the back side, offers clear and convenient controls. Easily adjust illumination levels for optimal clarity.



Mobile Base

With a mobile base that allows easy movement and positioning of the illumination light.

ABOUT US

Zhangjiagang Medi Medical Equipment Co.,Ltd, founded in 1995, with more than 20 years experience in developing hospital furniture, is one of the largest manufacturer of hospital bed, patient stretcher, medical cart and other related hospital furniture in China or even Asia.

Through scientific management complying with ISO9001&ISO13485, by adopting advanced facility like Germany TRUMPF CNC, Japan OTC Robot Welding, Auto-power coating line with Swiss GEMA gun , Medi focuses on providing brand quality product with CE approval.

Headquartered in Zhangjiagang nearby Shanghai, with subsidiary oversea offices in Hungary, Dubai, Thailand, New Zealand, Kenya, Colombia, Medi acting professionally and swiftly for both pre-sales and after-sales. The products are widely accepted by global clients, such as UK RIE hospital, UNIVERZITETNI KLINICNI CENTER, Kokilaben Dhirubhai Ambani Hospital, Singapore National University Hospital etc.



YOU TUBE



MEDIK

Pls scan the QR code for more info

www.mediksurgical.com | +86 139 6223 1292

Medik reserves the right to make changes without notice in design, specifications and models. The only warranty Medik makes is the express written warranty extended on the sale or rental of its products.

© 2022 Medik Services, Inc. ALL RIGHTS RESERVED. APR53703 rev 1 01-JAN-2022 ENG – US