

# YA-03F

Electric Field Surgery Table



*Better care  
Better comfort!*

# YA-03F

## ELECTRIC FIELD SURGERY TABLE

Minimized packing size and affordable carrying weight, necessary functions, all the features make the YA-03F electric operating table suitable for field operations like military or field rescue.

### FEATURES

- 1.Cold rolled steel with powder coating treatment
- 2.Lightweight but stable design
- 3.All operation via hand remote controller
- 4.Detachable design for small and quick loading
- 5.Suitable for field and special location application

ITEM	PARAMETER
Unfolded size	1980*470*(515~1050mm)
Folded size	960 * 550 * 530mm
Head board	+37 / -90 degree
Back board	-18 / +82 degree
Leg board	-31 / +56 degree
Lying surface	1980 * 470mm
Tren / Reverse Tren position	≥15 / 15 degree
Left and right tilt	≥15 / 15 degree
Height adjustment range	550 ~1050mm
Loading weight	≥135kg
Surgical table N.W	≤75kg

ITEM	PARAMETER
Mattress	√
Arm support	√
Leg support	√
Body support	√
Carrying box	√
Anesthesia screen	√
holder	√



**Back board adjustment**

Back board can be adjusted through hand remote controller from -18 to +82 degree.



**Head board adjustment**

Head board can be adjusted manually from -90 to +30 degree and removable.




**Tren position**

Tren position can be adjusted through hand remote controller from 0-15 degree.



**Reverse-Tren position**

Reverse-Tren position can be adjusted through hand remote controller from 0-15 degree.

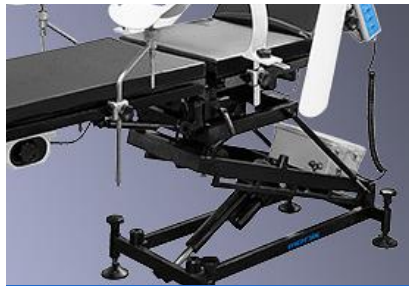

**Height adjustment**

Height can be adjusted through hand remote controller from 550 to 1050 mm.



**Leg board**

Leg board can be adjusted up and down from -31 to +56 degree by gas spring and outward adjustable manually, also is removable for gyne and obstetric application.



**Support holder**

Ground depth adjustment range can be adjusted via 4 pcs support holder.



**Left and right tilt**

Left and right tilt can be adjusted through hand remote controller from 0 to 15 degree.

## ABOUT US

Zhangjiagang Medi Medical Equipment Co.,Ltd, founded in 1995, with more than 20 years experience in developing hospital furniture, is one of the largest manufacturer of hospital bed, patient stretcher, medical cart and other related hospital furniture in China or even Asia.

Through scientific management complying with ISO9001&ISO13485, by adopting advanced facility like Germany TRUMPF CNC, Japan OTC Robot Welding, Auto-power coating line with Swiss GEMA gun , Medi focuses on providing brand quality product with CE approval.

Headquartered in Zhangjiagang nearby Shanghai, with subsidiary oversea offices in Hungary, Dubai, Thailand, New Zealand, Kenya, Colombia, Medi acting professionally and swiftly for both pre-sales and after-sales. The products are widely accepted by global clients, such as UK RIE hospital, UNIVERZITETNI KLINICNI CENTER, Kokilaben Dhirubhai Ambani Hospital, Singapore National University Hospital etc.



YOU TUBE



MEDIK

**Pls scan the QR code for more info**

**[www.mediksurgical.com](http://www.mediksurgical.com) | +86 139 6223 1292**

Medik reserves the right to make changes without notice in design, specifications and models. The only warranty Medik makes is the express written warranty extended on the sale or rental of its products.

© 2022 Medik Services, Inc. ALL RIGHTS RESERVED. APR53703 rev 1 01-JAN-2022 ENG – US