

YA-13EB

Carbon Fiber Spine Frame



*Better care
Better comfort!*

YA-13EB

CARBON FIBER SPINE FRAME

DESCRIPTION

YA-13EB mechanical spine frame is used for spinal surgery, cervical spine, lateral position and other surgical positions. It is suitable for any operating table. The all-carbon fiber material is convenient for C-arm and O-arm shooting. It can be adjusted in height, front and back, left and right. It can be disassembled and stored when not in use without taking up space. YA-13EB is the best choice for the operating room.

- FEATURES
- The spinal surgery frame is suitable for spinal surgery, cervical spine, lateral position, and other surgical positions.
 - The spinal surgery frame is suitable for any brand of operating table, through a simple connection device.
 - The full carbon fiber spinal support body is made of transparent material without any metal, which is convenient for C-arm and O-arm shooting without shadows.
 - The trend and reverse trend adjustment range of the spinal frame is ≥600mm, which can be used by children and adults.
 - The angle and position requirements of the spinal surgery frame can be adjusted through the ordinary operating table
 - The height of the spinal surgery frame can be adjusted by manual crank of the end column
 - The base uses four-point manual brake to ensure the safety and stability of the spine frame.
 - The spinal surgery frame can be disassembled and stored in a three-dimensional manner, which is convenient to place in any operating room or storage room without taking up space.


TECHNICAL
PARAMETER

Model	YA-13EB
Length	≥1400mm
Column Height	785-1195mm(±5mm)
Width	≥430mm
Patient Weight Capacity	300kg
Trend&reverse trend distance	600mm

ACCESSORIES

Name	Quantity
Main carbon fiber frame	1
Height adjustable base	1
Crank	1
Hand support	1
Angle support pad	4
X-ray translucent angle support	4
Carbon fiber clamp	2
Connecting adaptor	2
Net basket	1
Face support frame	1
Prone position head pad	1
Square support pad	1
Square support frame	1

DRAWING

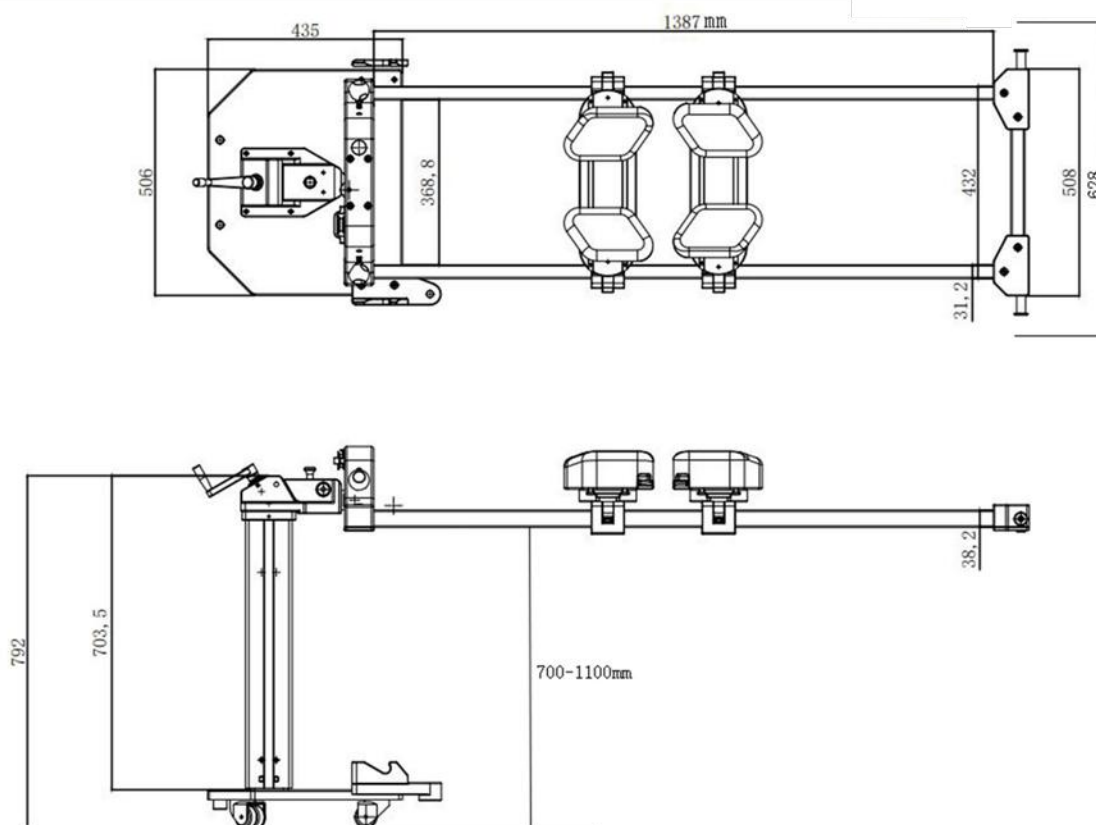



carbon fiber material

The full carbon fiber spinal support body is made of transparent material without any metal, which is convenient for C-arm and O-arm shooting without shadows.



Spinal surgery accessories

Usually comes with headrest, handrest and abdominal support pad.




Basket

The column is equipped with a basket to place items.



Mechanical crank and fixator

After unlocking the fixator, the height of the spinal frame can be adjusted by the crank to match the operating table in different positions.

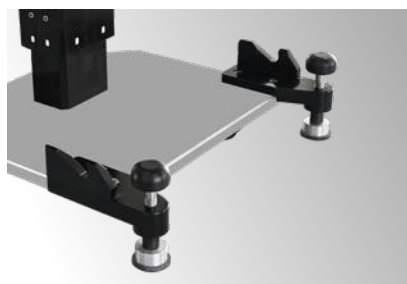

Brakes

The base uses four-point manual brake to ensure the safety and stability of the spine frame .



Angle indicator

The spine frame has left and right tilt angle indicator and trend angle indicator.

info@mediksurgical.com | +86 139 6223 1292

Medik reserves the right to make changes without notice in design, specifications and models. The only warranty Medik makes is the express written warranty extended on the sale or rental of its products.

© 2022 Medik Services, Inc. ALL RIGHTS RESERVED. APR53703 rev 1 01-JAN-2022 ENG – US