


YA-AT02

Autopsy Table

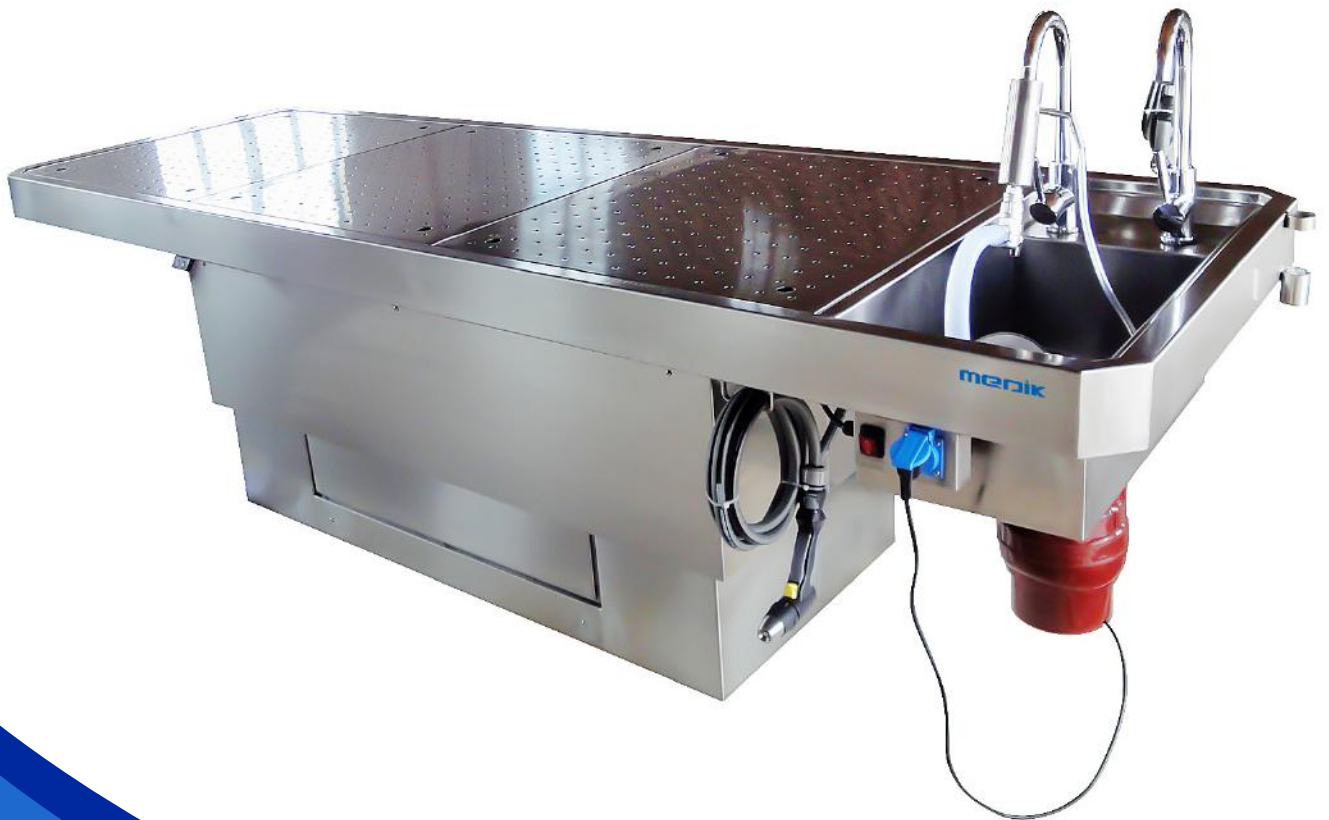

*Better care
Better comfort!*

YA-AT02

AUTOPSY TABLE

DESCRIPTION

Electrically height adjusting, special tap equipped for formalin absorbing, unique ventilation system for smell removal, whole SS304 constructed, all the superior features make the table as a preferred autopsy table for the teaching hospitals or medical universities.

FEATURES

- 1.High quality 304 stainless steel material, mirror polished surface
- 2.Double-layer structure, 3 sieved table-top plates.
- 3.Equipped with hot and cold water tap and liquid flushing faucets to effectively prevent virus cross-infection
- 4.The inside of the table is inclined downward to discharge harmful gases during the dissection process.
- 5.The crusher under the pool can crush broken tissues and small bones, and discharge them through the outlet after dissection.
- 6.With waterproof power sockets for convenient use of dissection tools
- 7.With a water spray gun with adjustable pressure to easily rinse away dirt.
- 8.Simple maintenance, safe and reliable, energy-saving, standard after-sales service, no worries.

STANDARD SPECIFICATION

Size	2600*860*750-1000mm
Height adjustment	750-1000mm
Accessory	water tap, crusher, power outlet, spray water gun
Capacity loading	350kg



Grid Plate

with 3 removable Grid Plates. Liquids and small objects may flow into the bottom through the holes.



Inclined Downward Table

The inside of the table is inclined downward to discharge harmful gases during dissection.



Water Tap

Equipped with 2 water taps, one tap can provide hot and cold water, and the other can flush formalin liquid.



Height Adjustment Button

2 buttons control the height adjustment of the autopsy table, with an overall stroke of 25cm.



Power Outlet

Equipped with waterproof power outlets for convenient use of dissection tools.



Crusher

There is a crusher under the pool that can be used to process and remove waste.



Spray Water Gun

Equipped with a water spray gun for cleaning dirt.

ABOUT US

Zhangjiagang Medi Medical Equipment Co.,Ltd, founded in 1995, with more than 20 years experience in developing hospital furniture, is one of the largest manufacturer of hospital bed, patient stretcher, medical cart and other related hospital furniture in China or even Asia.

Through scientific management complying with ISO9001&ISO13485, by adopting advanced facility like Germany TRUMPF CNC, Japan OTC Robot Welding, Auto-power coating line with Swiss GEMA gun , Medi focuses on providing brand quality product with CE approval.

Headquartered in Zhangjiagang nearby Shanghai, with subsidiary oversea offices in Hungary, Dubai, Thailand, New Zealand, Kenya, Colombia, Medi acting professionally and swiftly for both pre-sales and after-sales. The products are widely accepted by global clients, such as UK RIE hospital, UNIVERZITETNI KLINICNI CENTER, Kokilaben Dhirubhai Ambani Hospital, Singapore National University Hospital etc.



YOU TUBE



MEDIK

Pls scan the QR code for more info

www.mediksurgical.com | +86 139 6223 1292

Medik reserves the right to make changes without notice in design, specifications and models. The only warranty Medik makes is the express written warranty extended on the sale or rental of its products.

© 2022 Medik Services, Inc. ALL RIGHTS RESERVED. APR53703 rev 1 01-JAN-2022 ENG – US